


DEQAS October 2018 - 1,25-dihydroxyvitamin D Method Means * (\pm 1 SD)

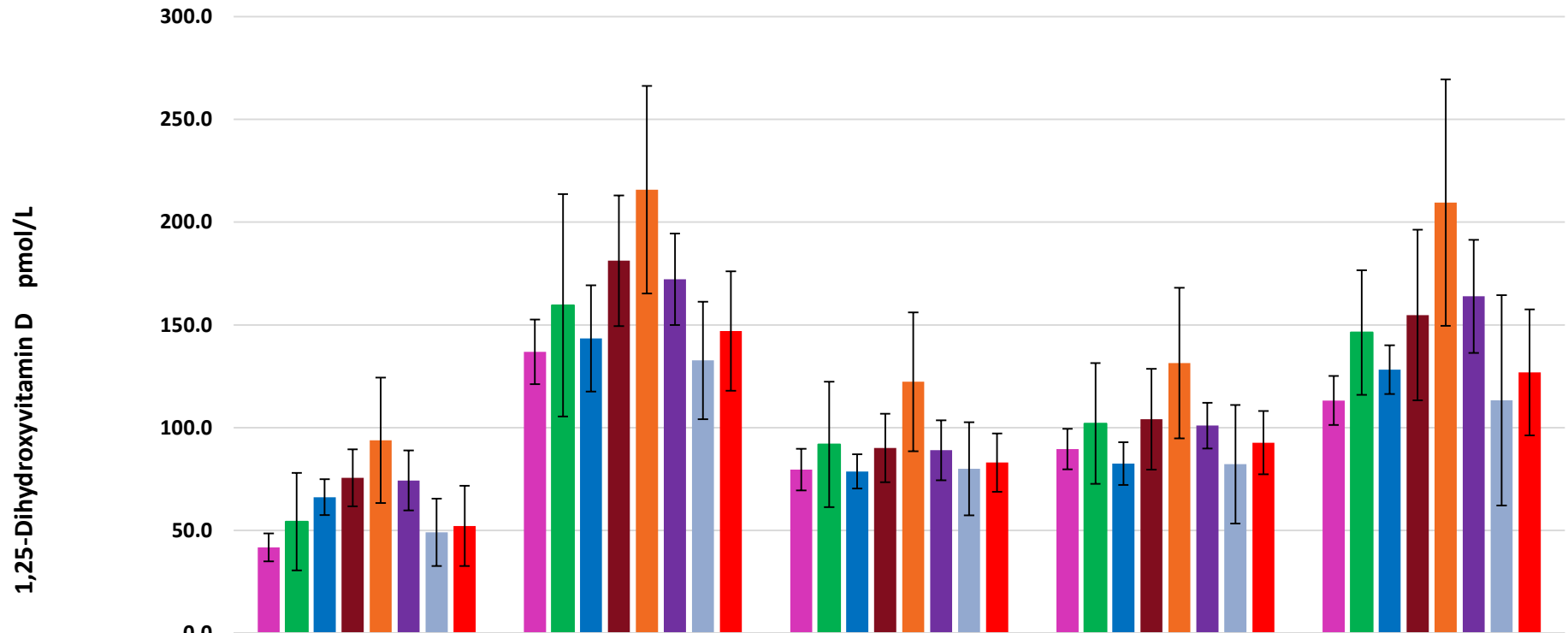

	Sample 431	Sample 432	Sample 433	Sample 434	Sample 435
■ DiaSorin Liaison XL (n=93)	41.7	136.9	79.6	89.6	113.2
■ DIASource CT (n=5)	54.3	159.5	91.9	102.0	146.3
■ IDS EIA (n=7)	66.2	143.4	78.7	82.5	128.2
■ IDS RIA (n=11)	75.6	181.2	90.1	104.1	154.8
■ IDS iSYS (n=10)	93.8	215.8	122.3	131.4	209.5
■ IDS iSYS New (n=19)	74.3	172.2	89.0	101.0	163.9
■ LC-MS/MS (n=13)	49.1	132.7	80.0	82.2	113.3
■ ALTM (n=158)	52.2	147.0	83.0	92.7	126.9

* Methods with 5 or more results returned