


DEQAS January 2018 - 25-OHD Method Means (+/-1SD) for Major Method Groups\*

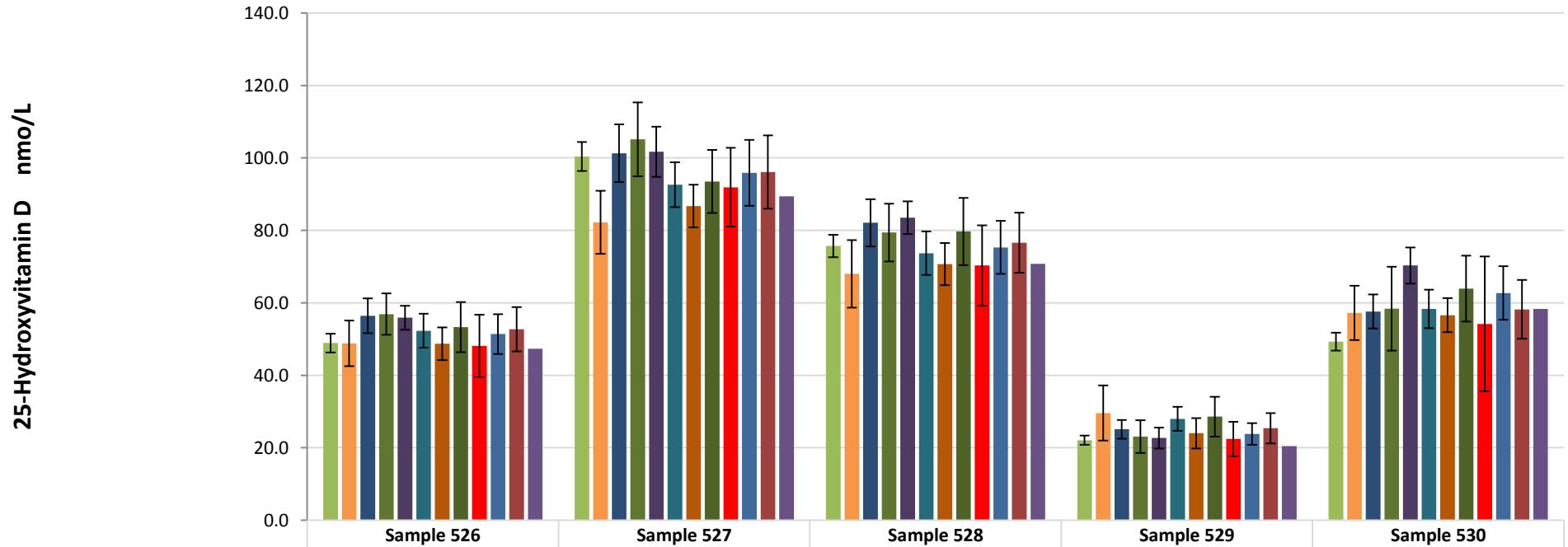

	Sample 526	Sample 527	Sample 528	Sample 529	Sample 530
Abbott Architect New Kit (n=70)	48.9	100.4	75.7	22.1	49.3
Beckman Unicel (n=32)	48.8	82.2	68.0	29.6	57.2
DiaSorin Liaison (n=197)	56.4	101.3	82.1	25.1	57.6
IDS iSYS (n=51)	56.9	105.1	79.4	23.1	58.4
IDS-iSYS New (n=14)	55.9	101.7	83.5	22.7	70.3
Roche Total (n=152)	52.3	92.6	73.7	28.0	58.3
Roche Vitamin D total II (n=25)	48.7	86.7	70.7	24.0	56.6
Siemens (n=65)	53.3	93.5	79.7	28.6	63.9
HPLC (n=16)	48.1	91.9	70.3	22.4	54.2
LC-MS/MS (n=154)	51.4	95.9	75.3	23.8	62.7
ALTM (n=842)	52.7	96.1	76.6	25.4	58.2
TARGET VALUE	47.3	89.4	70.8	20.5	58.3

\* Methods with 10 or more results returned