


DEQAS October 2017 - 25-OHD Method Means (+/-1SD) for Major Method Groups*

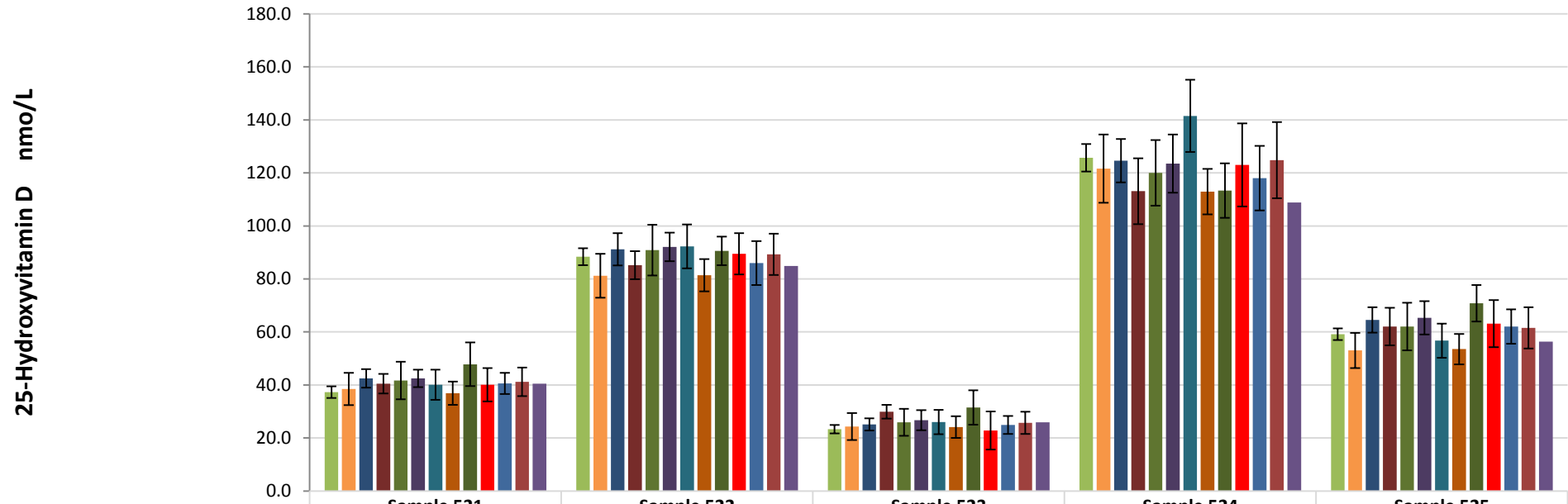

	Sample 521	Sample 522	Sample 523	Sample 524	Sample 525
Abbott Architect New Kit (n=70)	37.3	88.4	23.3	125.7	59.1
Beckman Unicel (n=34)	38.5	81.2	24.3	121.6	53.0
DiaSorin Liaison (n=209)	42.5	91.2	25.1	124.6	64.5
IDS EIA (n=12)	40.5	85.2	29.9	113.1	62.0
IDS iSYS (n=62)	41.7	90.9	25.9	120.0	62.0
IDS-iSYS New (n=10)	42.5	92.1	26.7	123.5	65.3
Roche Total (n=161)	40.1	92.3	26.0	141.5	56.7
Roche Vitamin D total II (n=18)	36.9	81.4	24.1	112.9	53.5
Siemens (n=64)	47.8	90.6	31.5	113.3	70.8
HPLC (n=19)	40.1	89.5	22.8	123.0	63.1
LC-MS/MS (n=159)	40.6	86.0	24.9	118.0	62.0
ALTM (n=871)	41.2	89.3	25.7	124.8	61.5
TARGET VALUE	40.5	84.9	25.9	108.8	56.3

* Methods with 10 or more results returned